

Miscellaneous Findings in Body Fluids

LEILANI COLLINS

ARDS = adult respiratory distress syndrome; BAL = bronchoalveolar lavage; LE = Lupus Erythematosus.

INDEX TERMS: body fluids; cerebrospinal; serous; synovial.

Clin Lab Sci 2009;22(1):54

LEARNING OBJECTIVES

1. Describe the significance of eosinophils in urine.
2. Describe lupus erythematosus cells in serous and synovial fluids.
3. Describe the significant cellular findings in bronchial alveolar lavage (BAL) specimens.

Leilani Collins MS MT(ASCP)SH CLS(NCA) is associate professor, Clinical Laboratory Science Program, University of Tennessee Health Sciences Center, Memphis TN.

Address for correspondence: Leilani Collins MS MT(ASCP)SH CLS(NCA), associate professor, Clinical Laboratory Science Program, University of Tennessee Health Sciences Center, 930 Madison Avenue, Suite 670, Memphis TN 38163. (901) 448-6299. lcollins@utmem.edu.

Leilani Collins MS MT(ASCP)SH CLS(NCA) is the Focus: Body Fluids guest editor.

Cytocentrifuge slides can provide important information in diagnosing various conditions.

Urine for eosinophils

Urine samples for eosinophils may be submitted when inter-

The Focus section seeks to publish relevant and timely continuing education for clinical laboratory practitioners. Section editors, topics, and authors are selected in advance to cover current areas of interest in each discipline. Readers can obtain continuing education credit (CE) through P.A.C.E.® by completing the continuing education registration form, recording answers to the examination, and mailing a photocopy of it with the appropriate fee to the address designated on the form. Suggestions for future Focus topics and authors, and manuscripts appropriate for CE credit are encouraged. Direct all inquiries to the Clin Lab Sci Editorial Office, IC Ink, 858 Saint Annes Drive, Iowa City IA 52245. (319) 354-3861, (319) 338-1016 (fax). ic.ink@mchsi.com

stitial nephritis secondary to antibiotic therapy is suspected. Slides to search for eosinophils should be prepared by concentration of the sample by either standard centrifugation or cytocentrifugation. Slides should be stained with Wright stain and a thorough search made for eosinophils. A finding of any eosinophils is significant.


Lupus Erythematosus (LE) cells

LE cells are neutrophils containing a smooth homogeneous nuclear mass. The nucleus of the neutrophil will be displaced to the edge of the cell wall by the large homogeneous mass. In patients who have systemic lupus erythematosus, the characteristic LE cell may be seen in serous and synovial fluids. All the factors necessary for formation of these cells—incubation, trauma to the cells, and the LE factor—are present *in vivo* in patients with systemic lupus erythematosus. The finding of these cells can be of diagnostic significance in patients not previously identified with this disease. See Figure 1.

Bronchoalveolar lavage (BAL)—the fluid that isn't a fluid

Bronchoalveolar lavage (BAL) is usually performed on hospitalized patients who have non-resolving pneumonia in an attempt to identify organisms that are not responding to antibiotic treatment or to detect malignancy. The “fluid” is obtained by introducing warmed saline into the lungs through

Figure 1. Two LE cells in synovial fluid

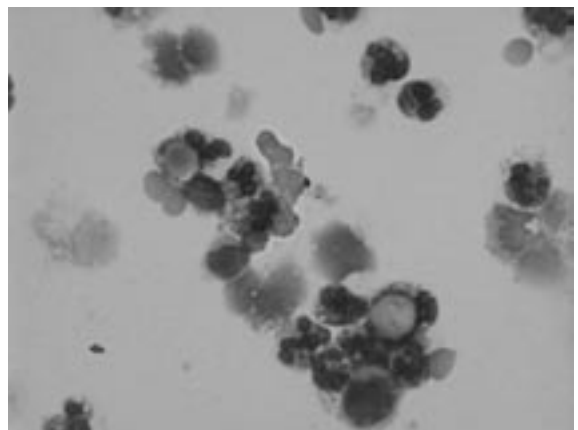

a bronchoscope in 20 mL aliquots followed by aspiration of the saline. A total of 100-120 mL of saline is introduced and aspirated with a recovery of 40-60 mL.³ The “fluid” is sent to the laboratory for extensive microbiology testing and cytology examination. Occasionally, cell counts and differential counts are requested. Any laboratory performing these procedures needs to be aware that organisms in these specimens are airborne and caution should be taken when containers of fluid are mixed and opened. Containers should be opened under a biological safety hood, tests should be “set up” under the hood (i.e., hemacytometers charged, cytofunnels filled, films prepared), and laboratory personnel performing these procedures should wear a mask while performing cell counts. There are no reference ranges for cell counts on these “fluids” and they have no clinical significance. As in any fluid, the types of cells present are significant. In BAL, the most clinically meaningful information is supplied by the extensive microbiological testing for bacteria and fungus that is the purpose for performing the procedure. The presence of numerous neutrophils can indicate infection or inflammation; numerous eosinophils can indicate allergic reactions, parasitic infections, or eosinophilic pneumonia; and excessive lymphocytes can indicate conditions including viral pneumonia, connective tissue disorders, or lymphoma.³ Macrophages containing tar can be seen in samples from patients who smoke. Erythrophages may be seen if there has been hemorrhage in the lungs more than 12 hours prior to the procedure. Ciliated epithelial cells are columnar epithelial cells found in the upper respiratory tract. They are elongated cells with the nucleus on one end of the cell and cilia on the other.


See Figure 2. If seen on the hemacytometer when performing cell counts on fresh specimens, these cells appear mobile as they are propelled by their cilia. Type II pneumocytes, resembling adenocarcinoma, can be seen in patients with adult respiratory distress syndrome (ARDS). Type II pneumocytes are very large, dark-staining cells in clumps containing large vacuoles that contain surfactant. Bacteria engulfed by neutrophils may be seen in cases of bacterial infection. Pneumocystis organisms, if present, will appear as amorphous material that, upon careful examination, may contain the round cysts characteristic of the organism. See Figure 3.

Clin Lab Sci encourages readers to respond with thoughts, questions, or comments regarding this Focus section. Email responses to george@fritsmafactor.com. In the subject line, please type “CLIN LAB SCI 22(1) FO BODY FLUIDS”. Selected responses will appear in the Dialogue and Discussion section in a future issue. Responses may be edited for length and clarity. We look forward to hearing from you.

WORKS CITED

1. Strasinger S, DiLorenzo M, editors. Urinalysis and body fluids, 5th edition, Philadelphia: F.A. Davis; 2008.
2. Collins L. Body fluids in the hematology laboratory. In: Rodak B, Fritsma G, Doig, K., editors. Hematology: clinical principles and applications, 3rd edition. St. Louis: Elsevier; 2007.
3. Mego T. Body fluids. Presented at CLSA Annual Meeting, Anchorage AK, April 2008. Available from <http://www.clsaonline.org/BODY%20FLUIDS%20Lecture%2004202008.pdf>. Accessed 2008 Sept 1.
4. Kjeldsberg C and Knight J, editors. Body fluids, 3rd edition. Chicago: ASCP Press; 1993.

Figure 2. Ciliated epithelial cells in BAL specimen

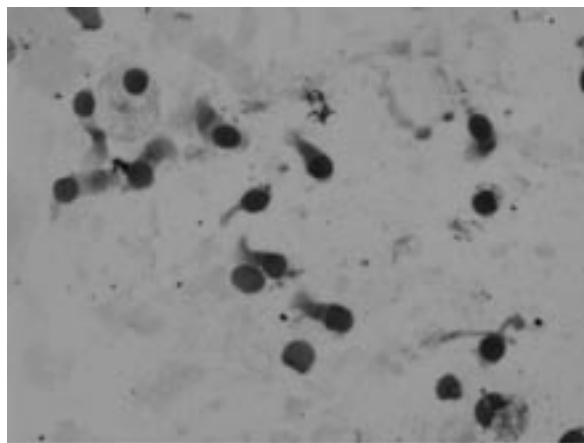

Note nucleus at one end of cell, elongated cytoplasm, and cilia at end opposite nucleus.

Figure 3. Pneumocystis in BAL specimen

